



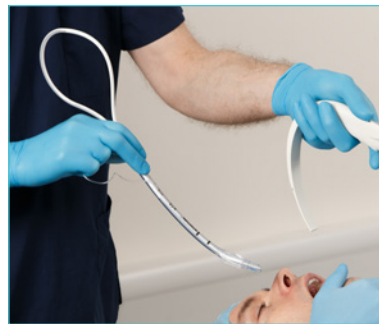

Universal Stylet Bougie

USB™

A hybrid stylet and bougie

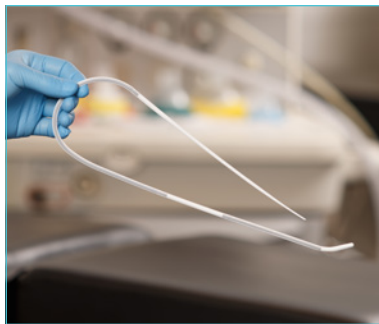

Airway Management ▪ Airway Accessories



Quality, innovation and trust

Interact with us



www.intersurgicalinc.com



Introducing the Universal Stylet Bougie (USB™)

A hybrid stylet and bougie in one simple to use device

The Universal Stylet Bougie, or USB, represents a significant advance in the design and development of tracheal tube introducers (bougies) and stylets.

Traditionally, the design requirements of stylets and bougies have always been considered incompatible. A stylet requires sufficient rigidity to allow the tracheal tube to be manipulated to and maintained at a variety of angles, whereas a bougie needs to be flexible enough to allow positive tracheal ring feedback while incorporating specific memory characteristics to allow correct placement.

The unique design of the USB incorporates two malleable sections with different and distinct properties on either side of a more flexible section. The USB can easily be manipulated to a variety of angles when used as a stylet, yet has the flexibility when used as a bougie to provide a significantly enhanced level of tracheal ring feedback.

In addition, the hexagonal shape provides less contact with the inner surface of the tracheal tube, providing particularly easy insertion and removal.

Key features


Low friction hexagonal shape

for optimal insertion and removal



Device information

confirmation of size, diameter and other key product details

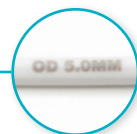

Flattened (non-patient) end

allows the user to insert the end of the USB back into the connector of the tracheal tube to retain product position

Coude tip

to optimize insertion. Atraumatic to reduce the potential for trauma

optiLube™

a sachet of lubricant is included in each pack to facilitate optimal insertion and removal (8090015 only)



Malleable sections and mechanical amplification

the sections of the device marked in grey indicate where the USB is malleable, allowing manipulation of the USB to the desired angle, and in combination with the flexible middle section, provide enhanced tracheal ring feedback

Length markings

confirmation of length at 45, 20 and 12 cm from the tip of the patient end of the device for correct positioning



Key benefits

The USB offers a number of key advantages over traditional bougies and stylets

Hybrid bougie and stylet — the unique design of the USB means the device can be used as a stylet or as a bougie.

Easier insertion and removal — the unique hexagonal shape of the USB provides less contact and interference with the tracheal tube ensuring the device can always be easily inserted and removed. The hexagonal shape also provides optimal grip and reduces the potential for rolling.

Easier storage — the folded design and smaller footprint of the USB makes it easier and more convenient to store than traditional stylets and bougies. It also ensures orientation of the device is always clear.

Improved tracheal ring detection — when there is a restricted laryngeal view, it is not always possible to visually determine whether a bougie has been correctly placed in to the trachea. In such circumstances, immediate feedback as to whether placement has been successful can be obtained by determining whether the tracheal rings can be felt as the bougie is moved backwards and forwards. This is the same for the USB, except its unique design, incorporating a flexible section positioned between two malleable sections, provides an enhanced level of amplification and tracheal ring detection not usually seen with traditional bougies.



Due to its unique design, the USB™ can be used in a variety of ways as both a bougie and a stylet

As a bougie:



As a standard bougie.



As a standard bougie, but with the USB retained in its original shape. The user can insert the device while still in this shape to retain directionality, then railroad the tube over the USB once correctly positioned.



The Tracheal Tube can be pre-loaded over the USB so the tip of the tube is located at the point where the second malleable section begins (identified by the grey sections). The proximal end is then folded back over into the Tracheal Tube so the tube does not move accidentally. Once the bougie has been correctly positioned, the Tracheal Tube is railroaded over the USB as normal.

As a stylet:



Straight to cuff.



Curved.

Using the USB™

Ideal for both pre-hospital and in-hospital use

The USB can be used in the pre-hospital as well as the in-hospital environment. In the pre-hospital setting, the folded design and smaller pack makes it easy to store in emergency vehicles and rescue backpacks or kits without any compromise to the integrity or shape of the device. The hybrid bougie and stylet design also eliminates the need to carry separate devices of each type, while retaining both clinical options.



Using the USB as a bougie

- Increased amplification, providing superior tactile ring feedback and detection.
- Hexagonal shape and folded design provides continuous confirmation of product orientation.
- Proximal end of the USB can be folded into the Tracheal Tube to prevent accidental movement until the point where the Tracheal Tube is to be railroaded over the USB.
- Can be used in a variety of configurations, allowing flexibility in use and the possibility to use the device without the need for an assistant.



Using the USB as a bougie

Using the USB as a stylet


- Hexagonal shape provides less contact with the Tracheal Tube and it is therefore easier to remove the USB from the Tracheal Tube once correct placement has been achieved.
- The flexible section between two malleable sections makes the overall device less rigid than many stylets, reducing the potential for trauma.



Using the USB as a stylet

Adjusting the USB to the desired angle

The USB has two malleable sections which can be bent to create the desired angle. The sections of the device marked in grey indicate where the USB is malleable.

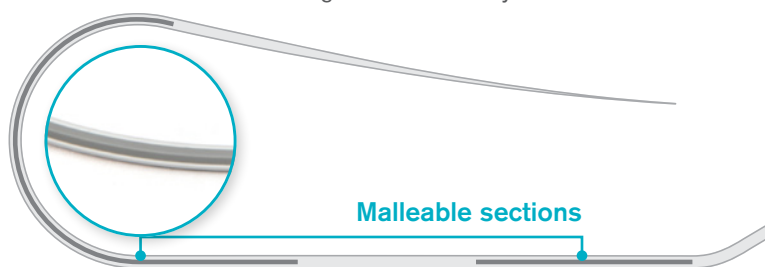

Code	Description	Qty.
8090015	Universal Stylet Bougie with OptiLube™	60 (Shipped in 60 units, 6 boxes of 10 units)
8080015	Universal Stylet Bougie	60 (Shipped in 60 units, 6 boxes of 10 units)

[Make an inquiry](#)
[Find out more](#)
[Watch the video](#)

OptiLube is a trademark of Optimum Medical.

U.S. Lit No. 2574 / IS6.16 US • Issue 7 09.25



Intersurgical Inc, 2055 Executive Drive, Indianapolis, IN 46241, USA
 T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • India • Ireland • Italy
 Japan • Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA