



## Sentri™ end tidal CO<sub>2</sub> mask

Delivers supplementary oxygen with CO<sub>2</sub> monitoring



**Oxygen and Aerosol Therapy ▪ Variable Oxygen Concentration (Low Flow)**



## Sentri™ Intersurgical EcoLite™ end tidal CO<sub>2</sub> mask

The Sentri adult mask permits the sampling of exhaled carbon dioxide in non-intubated patients during the administration of supplementary oxygen. There are also options available for connection to Microstream® capnography.

### Welcome to the Comfort Zone

The unique design and materials used to manufacture the Intersurgical EcoLite mask provides improved patient comfort and a reduced environmental impact.

### Comfortable for the Patient

Improved patient comfort has been key to the development of the Intersurgical EcoLite mask. The latest manufacturing technology has enabled us to combine two non-PVC materials in the same product, the polypropylene material forming the body of the Intersurgical EcoLite mask is clear, lightweight and rigid enough to maintain the masks shape, while the second softer TPE material is utilized in the seal, which is in contact with the patient's face.

### Comfortable for the Environment

The Intersurgical EcoLite mask is an important part of our Eco range, designed as part of our ongoing focus on sustainable development. Various groups have questioned the use of PVC in medical products and its impact on the environment. The utilization of alternative materials has resulted in the elimination of PVC from the Intersurgical EcoLite mask body and seal, reducing the environmental impact.

### Intersurgical EcoLite

[Find out more](#)

### Incurved nose seal

Conforms to different nose shapes designed to minimize the risk of oxygen entering patient's eyes

### Elastic can be positioned under or over the ears

Below ear position eliminates trauma to top of ears

### Soft face seals

Increased patient comfort

### Multi-channel oxygen tube

Oxygen still flows even if tube is kinked



### Luer lock port

Secure ETCO<sub>2</sub> monitoring line connection



### A choice of "under chin" or "on chin" positions

Provides a better fit on a wider range of patient face shapes

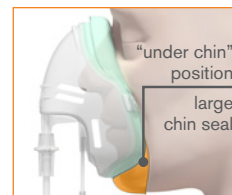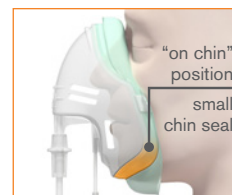

Code	Description	Tube length	Box Qty.	
1142015	Sentri Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line and tube	2.1 m	30	
1143015	Sentri Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line, filter and tube	2.1 m	30	
1143018	Sentri Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line, filter and tube with universal connector	2.1 m	30	
1141015	Sentri Intersurgical EcoLite, adult, mask with tube	2.1 m	30	
1142016	Sentri Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line for Microstream® capnography and tube	2.1 m	30	
1143017	Sentri Intersurgical EcoLite, adult, mask with CO <sub>2</sub> monitoring line for Microstream® capnography, filter and tube	2.1 m	30	

[Make an inquiry](#)

[Find out more](#)

Lower environmental impact product

Microstream is a registered trademarks of Oridion Medical 1987 Ltd.

**Please think before you print** Save energy and paper. If you must print this information sheet please print it double sided.

U.S. Lit No. 210.2 / IS10.24 US • Issue 17 09.24



Intersurgical Inc, 2055 Executive Dr. Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 [support@intersurgicalinc.com](mailto:support@intersurgicalinc.com) [www.intersurgicalinc.com](http://www.intersurgicalinc.com)

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • Ireland • Italy • Japan  
Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA