



## ClearLite™ anesthetic face masks

A clear vision for the future 



**Anesthesia** ▪ Anesthetic Face Masks



Quality, innovation and choice

Interact with us



[www.intersurgicalinc.com](http://www.intersurgicalinc.com)




## ClearLite™ anesthetic face masks

A range of anesthetic face masks which includes seven sizes to meet the clinical needs of all patient types and sizes. The range is available in sizes 0 neonatal to size 6 large adult.

Due to the enhanced anatomical design, the ClearLite mask offers superior seal performance with minimum pressure, reducing the risk of unnecessary discomfort or potential trauma.

The superior seal will also assist the clinician in providing optimum manual ventilation of the patient in both the anesthetic and emergency situation.

ClearLite is non-PVC and phthalate free which significantly minimizes the impact upon the environment. 

### Lightweight mask

Provides enhanced grip and comfort in use for the clinician

### Anatomical shape

Provides a precise fit for different face types and sizes. Minimum pressure required to achieve an effective seal reducing the risk of inadvertent trauma and nerve damage

### Crystal clear body

Provides excellent visibility for confirmation that the patient is well perfused

### Flexible non-inflatable cuff

Provides a comfortable and effective airtight seal

### Safe and secure connection


This is achieved between the mask and the breathing system in use. The ISO 22 mm tapered connection complies with ISO 5356

### Reduced mask volume

Reduces potential rebreathing and risks associated with elevated levels of end tidal CO<sub>2</sub>

### Color-coding

The comfortable, flexible seal is color coded to allow quick and easy identification of the correct sized mask

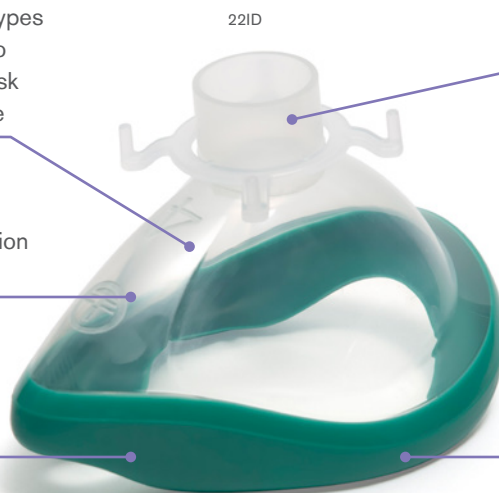

### Low environmental impact

Non-PVC and phthalate free

Code	Description	Color	Box Qty.
7296000	ClearLite, anesthetic face mask, size 6, extra large adult, red seal, 22ID		20 
7295000	ClearLite, anesthetic face mask, size 5, large adult, orange seal, 22ID		35 
7294000	ClearLite, anesthetic face mask, size 4, adult, green seal, 22ID		35 
7293000	ClearLite, anesthetic face mask, size 3, small adult, yellow seal, 22ID		35 
7292000	ClearLite, anesthetic face mask, size 2, paediatric, white seal, 22ID		25 
7291000	ClearLite, anesthetic face mask, size 1, infant, grey seal, no hook ring, 15OD		30 
7290000	ClearLite, anesthetic face mask, size 0, neonate, blue seal, no hook ring, 15OD		40 

[Make an inquiry](#)

[Find out more](#)

[Watch the video](#)

 Lower environmental impact product

**Please think before you print** Save energy and paper. If you must print this information sheet please print it double sided.

U.S. Lit No. 202.3 / IS5.31 US • Issue 4 09.24



Intersurgical Inc, 2055 Executive Drive, Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • Ireland • Italy • Japan  
Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA