




## Control-Cric™

Surgical cricothyrotomy kit

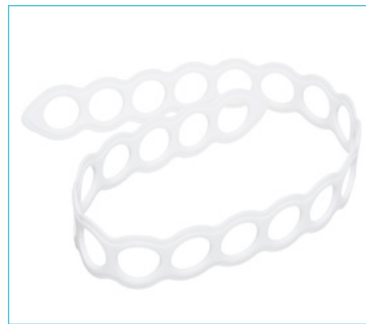

Airway Management ▪ Surgical cricothyrotomy



## Control-Cric™ surgical cricothyrotomy kit

The Control-Cric is packaged in 3 clear steps, guiding the provider through an otherwise chaotic procedure.

**Step 1** is the Cric-Knife™, a 10 mm dual-sided blade with an integrated sliding tracheal hook. This allows the provider to maintain control of the airway instead of reaching for additional tools and potentially losing the airway.

**Step 2** is the Cric-Key™, a 5.5 mm cuffed tube with an integrated introducer. The introducer gives tactile feedback of the tracheal rings, confirming placement.

**Step 3** is the neck strap, a flexible silicone strap to slide easily over the tube to maintain position during ventilation. A double swivel extension is also included with the kit to assist in minimizing tension and movement of the airway during patient care.



**Compact 3 step packaging**  
**Rapid deployment to oxygenate/ventilate**  
**BVM wedge included**



### Cric-Knife

- 10 mm dual-sided blade
- Integrated sliding tracheal hook
- Designed to grip like a pencil



### Cric-Key

- 5.5 mm cuffed tube
- Integrated introducer
- Provides tactile feedback of the tracheal rings to avoid false passage
- Airway tube is reinforced and is Radiopaque



### Neck strap

- Secure the tube after procedure
- Fits on 22 mm connection
- Wedge to assist with easy disconnect from the Cric-Key



### Extension tube

- Double swivel extension set
- Double flip top cap (Suction and Sealed Port to introduce a flexible scope)

| Code     | Description               | Box Qty.                                                                                |
|----------|---------------------------|-----------------------------------------------------------------------------------------|
| 313-4567 | Control-Cric kit sterile  | 5  |
| 313-4568 | Control-Cric training kit | 10                                                                                      |



Control-Cric, Cric-Key and Cric-Knife are manufactured by Pulmonary Inc and distributed by Intersurgical Ltd.

U.S. Lit No. 371.1 / IS6.21 US • Issue 1 04.26



Intersurgical Inc, 2055 Executive Drive, Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • India • Ireland • Italy  
Japan • Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA