



Cricalot

High-repetition task training kit for surgical cricothyrotomy



Airway Management ▪ Surgical cricothyrotomy



Cricalot

Cricalot is a high-repetition task trainer for surgical cricothyrotomy. Cricalot addresses the need for an anatomically accurate, realistic, low-cost training kit for cricothyrotomy, allowing clinicians to routinely practice this rarely occurring, potentially lifesaving procedure.

Feature and benefits

Cricalot has a range of features giving significant benefits to users.



Anatomically accurate

Each Cricalot kit contains a pair of models, one female and one male, with anatomically accurate thyroid cartilage, giving familiarity with both female and male anatomy. When opened, the outer case of each model acts as a sternum, allowing stabilization of the larynx and the operator's cutting hand, both critical elements of the cricothyrotomy procedure.

Skin analogue

200 skins and 2 rolls of cotton padding are provided with each Cricalot. Cotton padding can be placed into a skin pouch to achieve the desired thickness. Artificial blood may also be added to the skin pouch to create a more realistic model.

Easy to use and cost-effective

Skins are cost-effective and feature an adhesive back, allowing easy application to each model. After use, they are quick to remove and replace, with each model rapidly available to be used again.

Cricalot pack contents:


- 1 female model
- 1 male model
- 200 skins
- 2 rolls of cotton padding
- 6 spare tracheas

Sternum

The sliding outer case of each Cricalot model acts as a sternum when open, allowing stabilization of the larynx and operators cutting hand, both critical elements of the cricothyrotomy procedure

Skins

Easily changeable, cost-effective, plastic pouches allowing the addition of artificial blood for a more realistic model

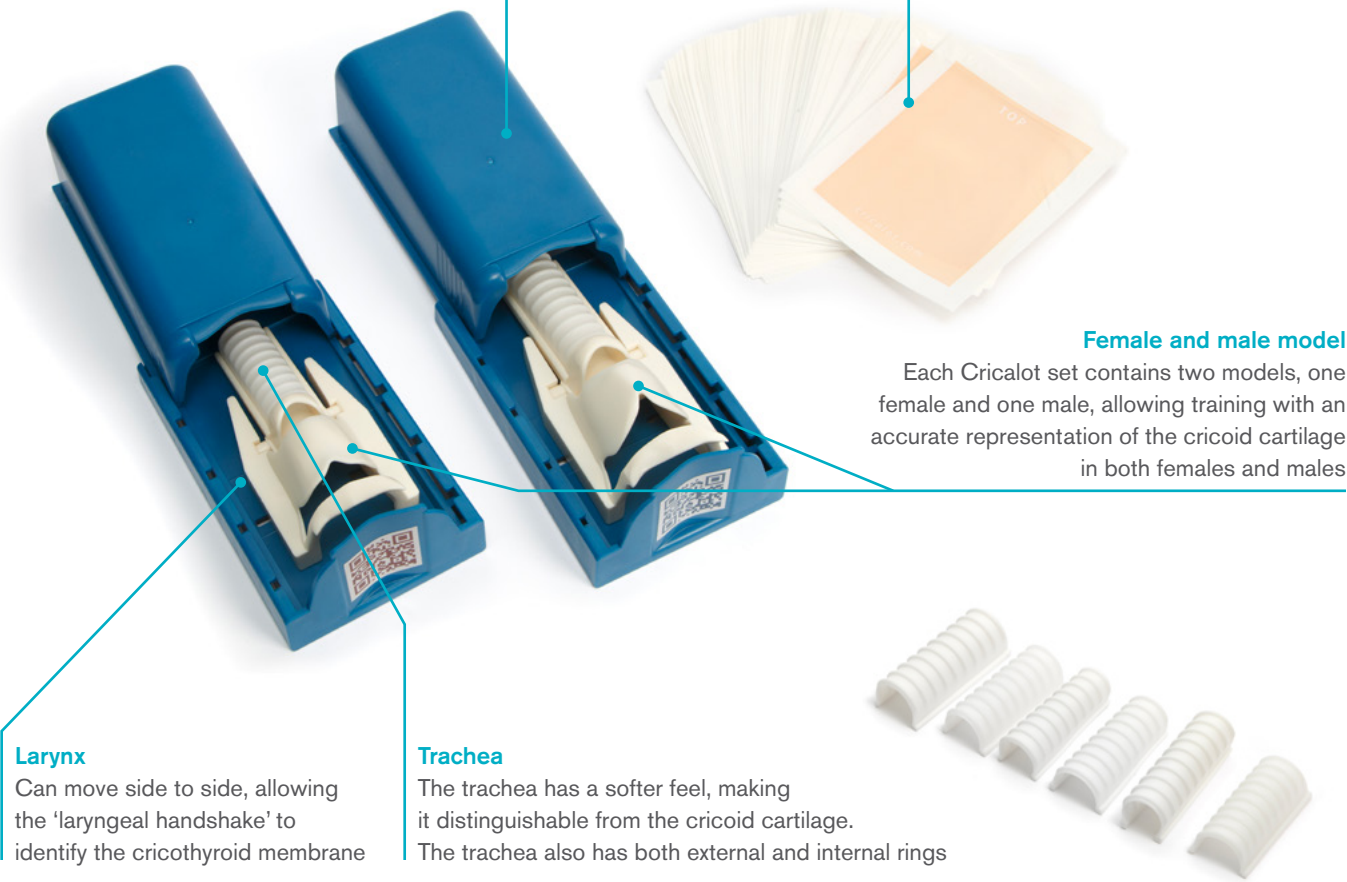

Code	Description	Box Qty.
3131030	Cricalot surgical cricothyrotomy trainer	5
3132030	Cricalot additional skins	1000



Cricalot is ideal for use with the Control-Cric™ training kit (Code: 313-4568)

The Control-Cric training kit is a non-sterile, training version of the Control-Cric, which allows familiarization with the Control-Cric and its features as part of a regular training program. The Control-Cric training kit ideally pairs with the Control-Cric with its distinctive outer packaging and yellow Cric-Key™ and Cric-Knife™.



Note: The Control-Cric training kit is not for human use.

Code	Description	Box Qty.
313-4568	Control-Cric training kit	10

Cricalot is manufactured by cricalot.com and distributed by Intersurgical Ltd. Control-Cric™, Cric-Key™ and Cric-Knife™ are manufactured by Pulmodyne Inc and distributed by Intersurgical Ltd.

U.S. Lit No. 372.1 / IS6.29 US • Issue 1 06.26



Intersurgical Inc, 2055 Executive Drive, Indianapolis, IN 46241, USA

T: 800-828-9633 F: 315-451-3696 support@intersurgicalinc.com www.intersurgicalinc.com

UK • Australia • Benelux • Brazil • Canada • China • Colombia • Czech Republic • Denmark • France • Germany • India • Ireland • Italy • Japan • Lithuania • Malaysia • Mexico • Norway • Philippines • Portugal • Russia • South Africa • Spain • Sweden • Taiwan • Turkey • USA